



The Acclaim 7000® — Multiple applications with a single laser

The Acclaim 7000 Nd:YAG laser allows you to expand your practice by offering cosmetic procedures such as:

- ▲ hair removal
- ▲ treatment of leg and facial veins
- ▲ our Dermique™ *Laser*FACIAL for sun-damaged skin

Advanced technology for hair removal.

The Acclaim 7000 offers selectable parameters for customizing treatment based on skin and hair type, improving outcomes on a variety of patients.

- ▲ Treats a wide variety of skin types (I-VI), including tanned skin.
- ▲ Provides high-peak power to effectively treat fine hair.

Treatment for a variety of vessels.

The Acclaim 7000 is the treatment of choice for a variety of vessels.

- ▲ Effective treatment of facial and leg telangiectasias (up to 3 mm).
- ▲ Successfully treats a wide variety of vascular lesions.
- ▲ Selectable laser parameters and spot sizes allow for a targeted approach to all veins.

Dermique *Laser*FACIAL.

The Acclaim 7000 Dermique *Laser*FACIAL promotes healthier looking skin and improves the appearance of sun-damaged skin.

- ▲ Reduces diffuse redness (the blush of rosacea) and minor facial veins.
- ▲ Quick and easy procedure with no patient downtime.

Increased productivity and profitability for your practice.

The Acclaim 7000 offers multiple applications, faster treatments, and greater patient throughput, making your practice more productive and profitable.

The Cynosure Acclaim 7000 laser.


*Another reason why we're
light-years ahead of the competition.*



“ The Acclaim 7000 has become vital to our practice. We not only use the laser for hair removal on all types of skin, but also for the treatment of facial and leg veins. Our patients are very satisfied with the results. ”

Susan Touma, MD

Before

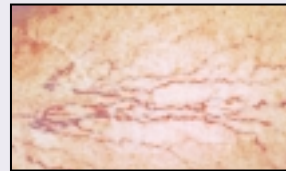

Leg vein removal

After



3 Mo Post 1 Tx




Facial vein removal



3 Mo Post 1 Tx


Courtesy Emil Tanghetti, MD

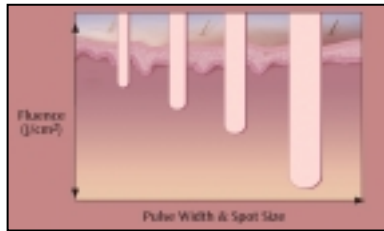

Short pulse durations and smaller spot sizes effectively treat smaller vessels. Deeper, larger vessels require longer pulse width and larger spot sizes.



The Acclaim 7000 offers a wide range of pulse widths to treat a broad variety of skin and hair types.

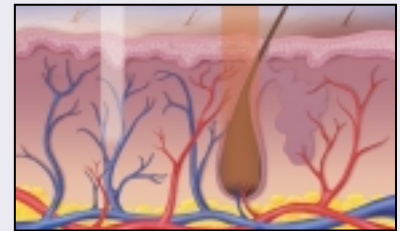

Effective treatment for hair and vein removal.

“ I have 20 lasers and light sources in my practice, but I did not realize how vital the Acclaim 7000 would be to my patients. We are achieving excellent results removing unwanted hair as well as treating facial and leg veins in patients who were not being adequately or safely treated with other devices. ”

Emil Tanghetti, MD



Facial hair removal




6 Mo Post 6 Tx

Courtesy Thomas Woo, MD



Leg hair removal

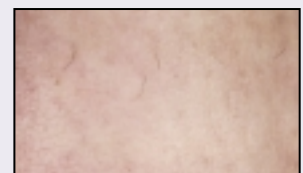

8 Mo Post 7 Tx

Courtesy of S. Galiberti MD



Leg vein removal



4 Mo Post 1 Tx

Corporate Headquarters

Cynosure, Inc.
10 Elizabeth Drive
Chelmsford, MA 01824
USA
Tel: 978.256.4200
Toll-Free: 800.886.2966
Fax: 978.256.6556

International Offices

People's Republic of China

Suzhou Cynosure Medical Devices, Co., Ltd.
5th Floor, Yuan Dong Da Sha
251 Chang Xu Road
Suzhou 215008, Jiangsu
People's Republic of China
Tel: +86 (512) 655 17051
Tel: +86 (512) 655 17351
Fax: +86 (512) 655 78485

Singapore

Cynosure Asia Pacific
220 Tagore Lane
#01-01/2 Liberty Warehouse
Singapore 787600
Tel: (65) 9688-4565
Fax: (65) 6853-4309

United Kingdom

Cynosure UK Ltd.
The Farm Offices
Lower Mount Farm, Long Lane
Cookham, Berkshire SL6 9EE
England
Tel: +44-1628 522252
Fax: +44-1628 520525

France

Cynosure SARL, Ltd.
86, Avenue Lenine
94250 Gentilly
France
Tel: +33-1-49.85.6005
Fax: +33-1-49.85.6004

Germany

Cynosure GmbH
Robert-Bosch-Str. 11-11B
D-63225 Langen, Germany
Tel: +49 6103 2011100
Fax: +49 6103 2011111

Japan

Cynosure KK
25-6 Taito 2-chome, Taito-Ku
Tokyo 150-0031, Japan
Tel: +81-3-5807-2761
Fax: +81-3-5807-2762

800-886-2966
www.cynosurelaser.com

PRODUCT SPECIFICATIONS

Feature	Acclaim 7000™
Wavelength (nm)	1064
Pulse Width (msec)	adjustable 0.4 - 300 —
Maximum Fluence per Spot Size	12 ² 50 J/cm ² 10 80 J/cm ²
Independent of pulse width, except 0.4 & 1 msec	7 160 J/cm ² 5 240 J/cm ² 3 300 J/cm ²
Maximum Rep Rate Adjustable, Continuous Pulsing or Single Pulse Operation	5Hz
Aiming Beam	Blue/Green
Skin Cooling System	Cold Air or Integrated Cooling
Electrical Requirements	110 VAC/20A ¹ or 220 VAC/30A
Size (HxWxD)	41" x 15" x 25" (104cm x 38cm x 64cm)
Weight	180 lbs (82 kg)

¹ reduced repetition rate

² optional accessory



 **CYNOSURE**[®]
Where art and science meet.