



The Apogee Elite™ Two-in-One Laser System — Innovative technology from Cynosure

Now you don't have to choose. The Apogee Elite contains two laser systems – a 755 nm Alexandrite laser and a high powered 1064 nm Nd:YAG laser – in one convenient unit. The Elite offers the most optimal wavelengths for:

- ▲ hair removal
- ▲ pigmented lesions
- ▲ facial and leg telangiectasias
- ▲ treatment of sun-damaged skin

Twice as effective, twice as easy to use.

- ▲ **Hair removal** for all skin types. 755 nm Alexandrite for skin types I-IV; 1064 nm Nd:YAG for skin types V-VI and tanned skin.
- ▲ **Pigmented lesions.** The Apogee Elite's short pulse durations provide effective treatment of pigmented lesions, age spots, freckles and sun spots.
- ▲ **Facial and leg telangiectasias.** Effective treatment for vessels up to 3 mm.
- ▲ **Facial treatment for sun-damaged skin.** Combination treatment visibly reduces redness and pigmentation, resulting in healthier, more radiant looking skin.

Apogee Elite advantages:

Efficacy, selectivity and speed.

- ▲ Choose the wavelength that provides optimal results for the removal of unwanted hair.
- ▲ Pulse duration selectivity allows you to customize treatments according to your patients' aesthetic needs.
- ▲ The Apogee Elite delivers high fluences from large spot sizes ensuring faster, more effective treatment.

Increased productivity and profitability for your practice.

The Apogee Elite laser system offers multiple applications and greater patient throughput, making your practice more productive and more profitable.

The Cynosure Apogee Elite laser system.

Another reason why we're light-years ahead of the competition.



“ I have used the Cynosure Alexandrite lasers since 1997. I have found them to be the most effective, safest and fastest lasers for hair removal. My patients have been consistently satisfied with their results while having minimal, if any, downtime associated with their treatments. I have been able to treat all skin types with this laser but have found the Nd:YAG to have a significant advantage with the darker skin types IV-VI as well as tanned skin. In addition, the Nd:YAG has been unmatched in ease-of-use and effectiveness in the treatment of both facial and lower extremity telangiectasias. ”

Peter Klainer, MD

Before



Facial hair removal

After



Courtesy Brooke A. Jackson, MD

9 Mo Post 7 Tx



Back hair removal




Courtesy Christopher Nanni, MD

7 Mo Post 4 Tx



Pigmented lesions

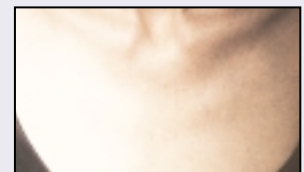


Courtesy Christopher Nanni, MD

2 Mo Post 1 Tx



Facial hair removal

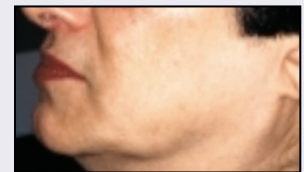


Courtesy Brooke A. Jackson, MD

8 Mo Post 6 Tx



The Apogee Elite offers a wide range of pulse widths to treat a variety of skin and hair types.

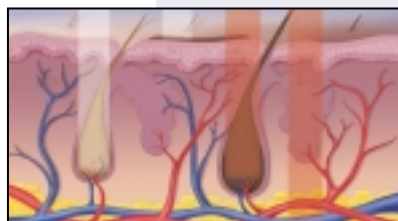

Choice of wavelengths effectively treats unwanted hair, veins, pigmentation, and improves skin tone.



Short pulse duration and small spot sizes effectively treat smaller vessels. Deeper, larger vessels require longer pulse widths and larger spot sizes.



Facial vein removal



Courtesy Emil Tanghetti, MD

2 Mo Post 2 Tx

Corporate Headquarters

Cynosure, Inc.
10 Elizabeth Drive
Chelmsford, MA 01824
USA
Tel: 978.256.4200
Toll-Free: 800.886.2966
Fax: 978.256.6556

International Offices

People's Republic of China

Suzhou Cynosure Medical Devices, Co., Ltd.
5th Floor, Yuan Dong Da Sha
251 Chang Xu Road
Suzhou 215008, Jiangsu
People's Republic of China
Tel: +86 (512) 655 17051
Tel: +86 (512) 655 17351
Fax: +86 (512) 655 78485

Singapore

Cynosure Asia Pacific
220 Tagore Lane
#01-01/2 Liberty Warehouse
Singapore 787600
Tel: (65) 9688-4565
Fax: (65) 6853-4309

United Kingdom

Cynosure UK Ltd.
The Farm Offices
Lower Mount Farm, Long Lane
Cookham, Berkshire SL6 9EE
England
Tel: +44-1628 522252
Fax: +44-1628 520525

France

Cynosure SARL, Ltd.
86, Avenue Lenine
94250 Gentilly
France
Tel: +33-1-49.85.6005
Fax: +33-1-49.85.6004

Germany

Cynosure GmbH
Robert-Bosch-Str. 11-11B
D-63225 Langen, Germany
Tel: +49 6103 2011100
Fax: +49 6103 2011111

Japan

Cynosure KK
25-6 Taito 2-chome, Taito-Ku
Tokyo 150-0031, Japan
Tel: +81-3-5807-2761
Fax: +81-3-5807-2762

PRODUCT SPECIFICATIONS

Feature	Apogee Elite™	
Wavelength (nm)	755	1064
Pulse Width (msec)	adjustable 0.5 - 300 —	adjustable 0.4 - 300 —
Maximum Fluence per Spot Size Independent of pulse width	15	25 J/cm ²
	12	35 J/cm ²
	10	50 J/cm ²
	7	100 J/cm ²
	5 ¹	25 J/cm ²
	3	—
300 J/cm ²		
Maximum Rep Rate Adjustable, Continuous Pulsing or Single Pulse Operation	2Hz	5Hz
Aiming Beam	Blue/Green	
Skin Cooling System	Cold Air or Integrated Cooling	
Electrical Requirements	110 VAC/20A ² or 200/220 VAC/30A	
Size (HxWxD)	41" x 15" x 25" (104cm x 38cm x 64cm)	
Weight	180 lbs (82 kg)	

¹ available with 0.5msec pulse width only, maximum rep rate 5Hz on the Alexandrite laser

² reduced repetition rate



800-886-2966
www.cynosurelaser.com

 **CYNOSURE**[®]
Where art and science meet.